



Tapestri Platform for Single-cell Multi-omics

Uncover genotype and phenotype insights simultaneously from single cells

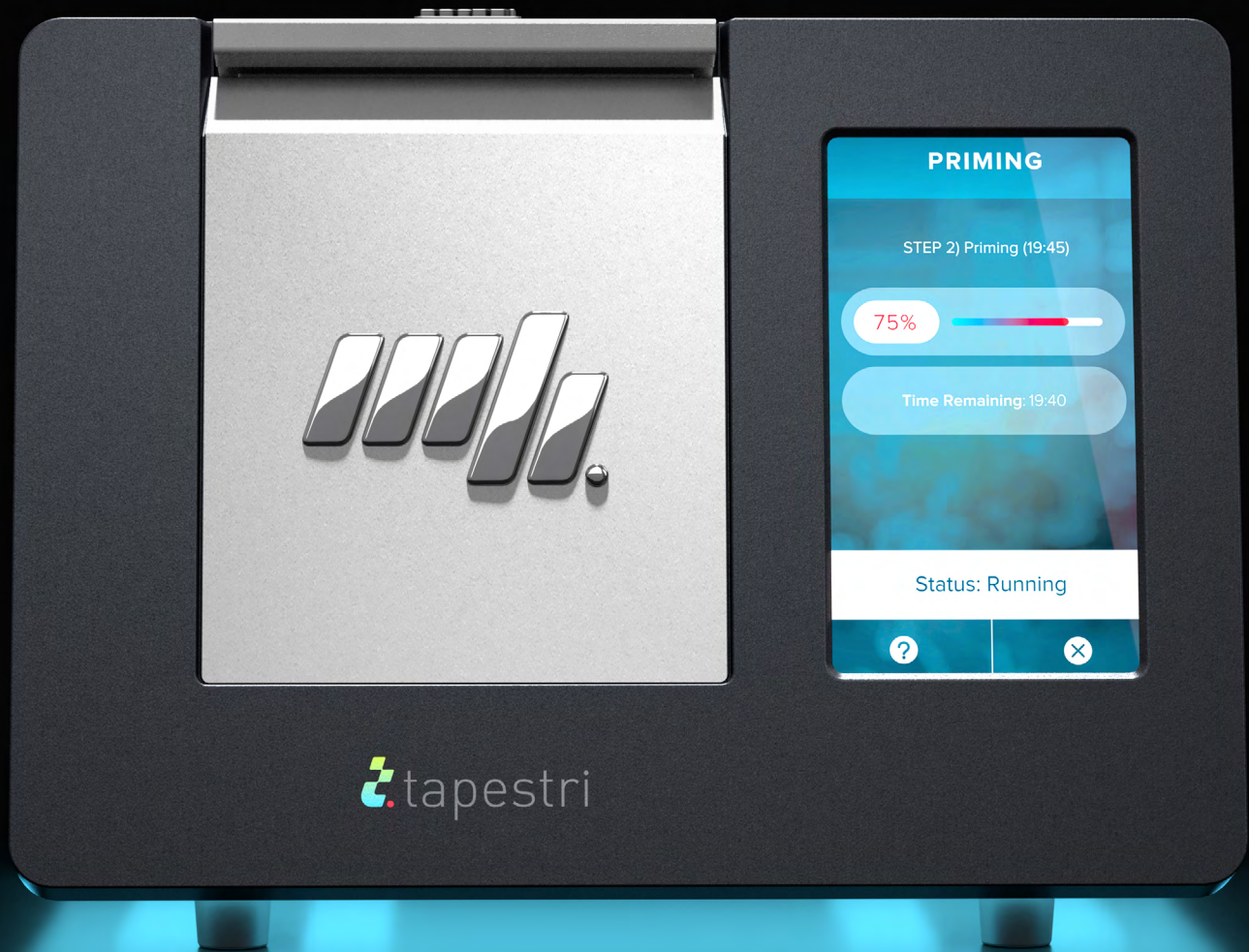

Unravel Cancer's Complexity with True Single-cell Multi-omics

The Tapestri® Platform is the world's first and only single-cell solution that provides both genotype and phenotype data from the same cell. Based on a novel two-step microfluidic workflow, the Tapestri Platform can access DNA and protein in single cells to give you a true multi-omics picture. With its unprecedented speed and scale, you can now get actionable insights needed to get to the ground truth and unravel the threads of cancer's complexity.

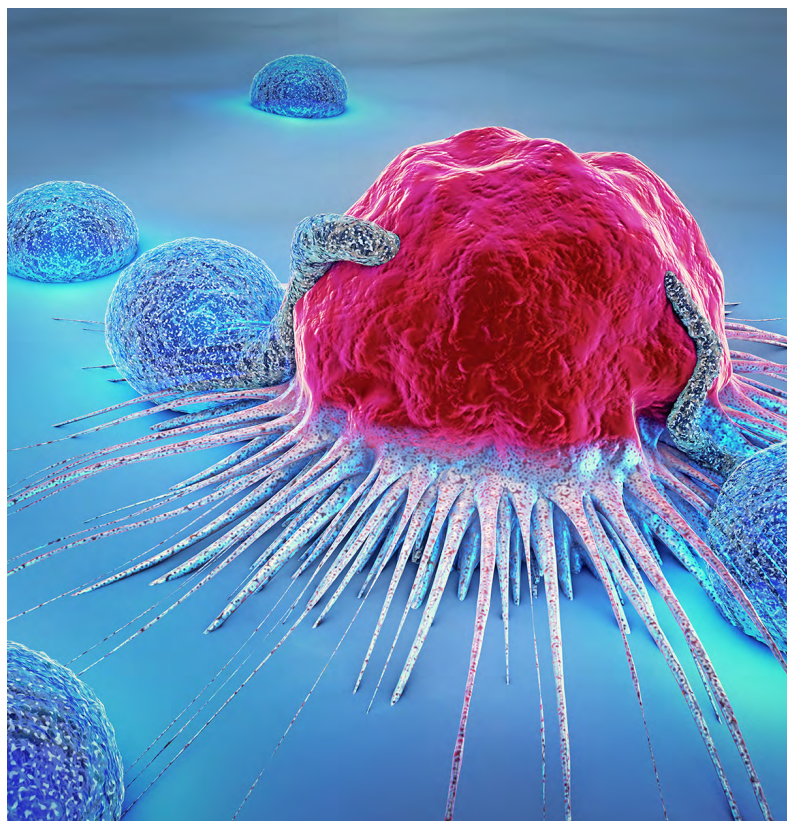


Tapestri Platform Highlights



Analyze DNA and protein expression simultaneously in single cells for true multi-omic insight



Complete solution with core kits to take you from single cells to a sequencing-ready library and analysis tools to convert your multi-analyte data to actionable insights

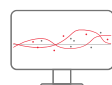

Targeted, customizable content for key oncology applications


The Importance of Understanding Cancer at a Single-cell Level

Cancer is a heterogeneous disease. In order to resolve heterogeneity and improve patient stratification, therapy selection, and disease monitoring, you need a tool that provides a holistic view of a cell, integrating data across multiple dimensions.

Unlike bulk sequencing, single-cell multi-omics enables you to:


- Detect rare cell populations
- Measure zygosity
- Identify co-occurring mutations
- Resolve clonal heterogeneity

CLONAL RESOLUTION WITH SINGLE-CELL PRECISION

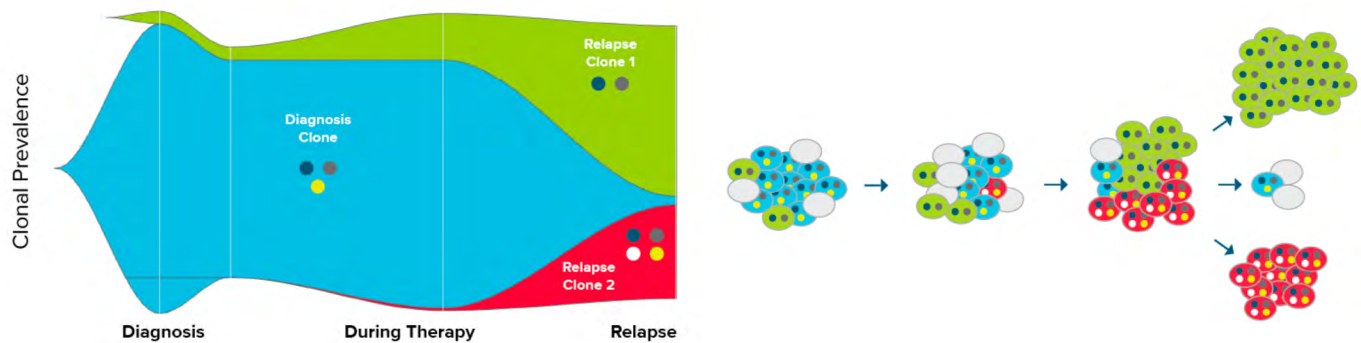

In order to see the evolution you need to understand the disease at a single-cell level.

THE POWER OF SINGLE-CELL MULTI-OMICS

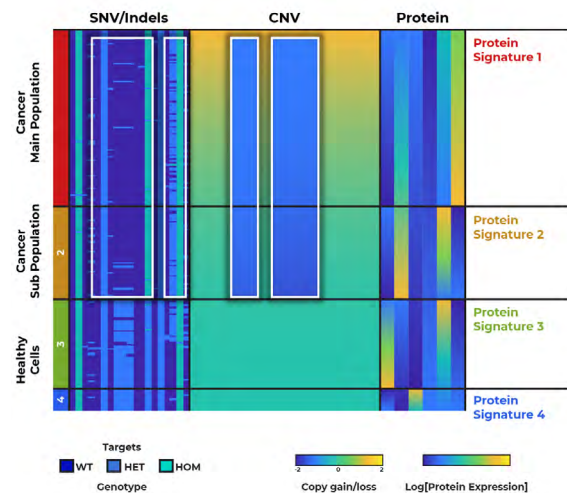

Combining genotype and phenotype data from single cells offers the resolution for uncovering unique disease signatures for personalized therapeutics.

“A single cell multiomics platform really gives you a degree of detail, which we have inferred or assumed, but have never been able to prove until now. What we can do is drill down to confirm that those [DNA mutations and immunophenotypes] co-exist . . . what is going to make a difference for that patient, it’s being able to target what is the true driving cause of the disease.”






- Katherine D. Cummins, MD, FRACP, FRCPA



Targeted, Customizable Content for Key Applications

Tapestri® Single-cell DNA Panels allow you to focus on the mutations and regions of interest that are most relevant to your disease research. To get your experiments running with minimal time and effort, consider our pre-designed and wet lab-tested Tapestri Single-cell DNA Catalog Panels. Alternatively, choose from our new Tapestri Single-cell DNA Published Panels designed by leading researchers and verified for performance through peer-reviewed publications. If you require more flexibility, consider our pre-designed but customizable Tapestri Single-cell DNA Virtual Design Panels. For maximum versatility, create a Tapestri Single-cell DNA Custom Panel. To investigate genotype and phenotype across thousands of individual cells, add a Protein Panel (oligo-conjugated antibodies from BioLegend) into your Tapestri experiments.


Applications include hematologic malignancies, solid tumors, genome editing, biomarker discovery and cell and gene therapy.

 Hematologic malignancies	 Solid tumors	 Genome editing	 Biomarker discovery	 Cell & gene therapy
<p>Catalog Panels:</p> <ul style="list-style-type: none"> Acute myeloid leukemia Chronic Lymphocytic Leukemia Myeloid <p>Published Panels:</p> <ul style="list-style-type: none"> Acute Lymphoblastic Leukemia Myeloid Clonal Evolution Myeloproliferative Neoplasms <p>Virtual Design Panels:</p> <ul style="list-style-type: none"> Acute Lymphoblastic Leukemia Chronic Myeloid Leukemia Classic Hodgkin's Lymphoma Diffuse Large B-Cell Lymphoma Follicular Lymphoma Mantle Cell Lymphoma Multiple Myeloma Myelodysplastic Syndromes Myeloproliferative Neoplasms T-Cell Lymphoma <p>Protein Panel:</p> <ul style="list-style-type: none"> TotalSeq-D Heme Oncology Cocktail 	<p>Catalog Panel:</p> <ul style="list-style-type: none"> Tumor hotspot 	<p>Partner with Mission Bio to implement single-cell multi-omics capabilities for genome editing, biomarker discovery and cell & gene therapy applications, and build a single-cell DNA or protein panel for the Tapestri Platform.</p> <p>Contact our Pharma Assay Development (PAD) representative for more information.</p>		
<p>Custom DNA Panels: 20 to 1,000 amplicons to cover DNA regions of interest</p>		<p>Custom Protein Panels: Configure an antibody cocktail using TotalSeq-D Catalog Antibody Oligonucleotide Conjugates</p>		

End-to-end Solution that Seamlessly Plugs into your NGS Workflow

Use the Tapestri instrument, reagents and consumables up-front of your NGS system and then Tapestri Pipeline and Tapestri Insights software for data analysis and visualization.


THE TAPESTRI WORKFLOW

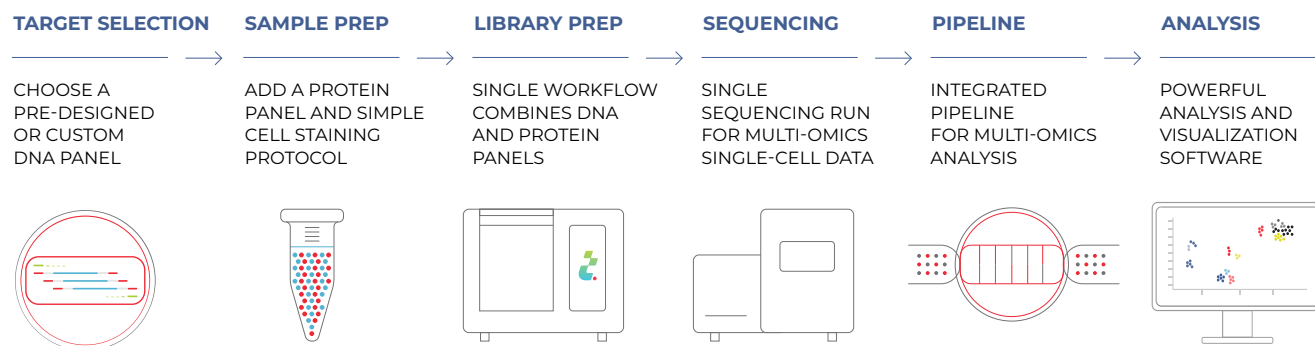

From Complex Multi-analyte Data to Actionable, Ground-truth Insights

Tapestri Pipeline and Tapestri Insights software solutions provide a streamlined bioinformatics workflow optimized for single-cell DNA and protein analysis. From sequence import to data analysis and visualization all packaged in a user-friendly experience, our turnkey analysis solutions ensure that you gain meaningful insights to advance your research.

TWO STEP PROCESS

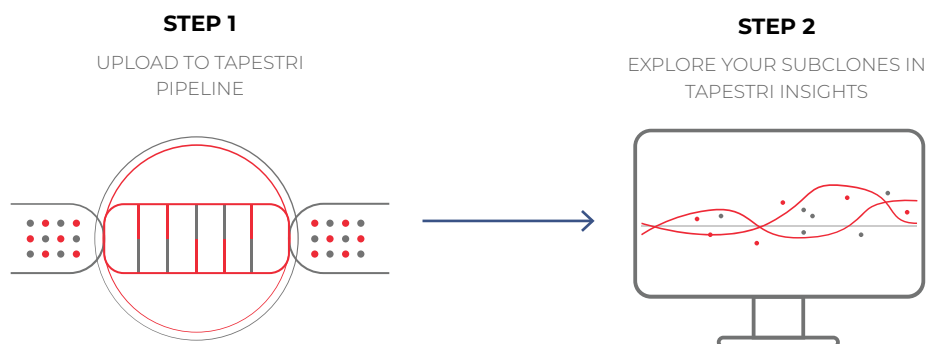

Tapestri Platform Specifications

TAPESTRI PLATFORM SPECIFICATIONS

Variant and analyte type detected	SNV, indel, CNV, LOH and translocation in DNA; protein expression
Cell input requirement	20,000-100,000 cells
Throughput	Up to 10,000 cells
Instrument dimensions	In H x W x D, 12.5 in x 11.75 in x 12.25 in (31.75 cm x 29.85 cm x 31.33 cm)
Operating conditions	Standard laboratory environments. Not exceeding 6,562 ft (2,000 m) above sea level. Maintain 5 – 85 % relative humidity, non-condensing. Maintain 4 in (10.2 cm) of clearance at the vents (back of instrument).
Power requirements	100 – 240 V AC power at 50 – 60 Hz; 2.0 A maximum

TAPESTRI PLATFORM COMPONENTS PART NUMBER

Tapestri Instrument	191335
Tapestri Single-Cell DNA Cartridge Kit	046459
Tapestri Single-Cell DNA + Protein Core Kit	MB03-0034
Tapestri Single-Cell DNA Panel Kits	missionbio.com/panels
Tapestri Single-Cell DNA Custom Panel Kits	missionbio.com/panels/custom-panels/

CONTACT US

Mission Bio, Inc.
6000 Shoreline Court, Suite 104
South San Francisco, CA 94080
+ 1.415.854.0058

info@missionbio.com
www.missionbio.com



For Research Use Only. Not for use in diagnostic procedures.
© 2021 Mission Bio, Inc. All rights reserved.
Brochure_MissionBio_Tapestri_Platform_RevG